




Photostability Testing for Cosmetics & Pharmaceuticals (ICH Guidelines & Forced Degradation)


The Q-SUN® xenon test chamber from Q-Lab offers state of the art, full-spectrum photostability test exposures for pharmaceutical and cosmetics applications. Q-SUN testers are available in two sizes with options that meet ICH Guidelines for Photostability Testing of New Drug Substances and Products, including simulation of direct sunlight or sunlight through window glass. Accurate control of light intensity and temperature ensures repeatable and reproducible results. These features enable the Q-SUN chamber to perform photostability tests according to in-vitro test procedures for sunscreens.

Key Features: Simplicity & Accuracy

- Simple operation and lamp replacement so you can concentrate on what's most important: your tests!
- Full spectrum xenon arc lamp and optical filter system creates D65 and ID65 sunlight simulations that enable you to test your products against outdoor or indoor daylight.
- SOLAR EYE® Irradiance Control monitors and controls the light intensity to ensure your tests are repeatable and reproducible.
- Precise control of chamber air or black panel temperature gives you flexibility to determine the effects of different temperature environments on photodegradation.
- Optional chiller unit enables tests to be run at room temperature or as low as 15°C so you can determine the synergistic effects of light and heat on your materials and products.

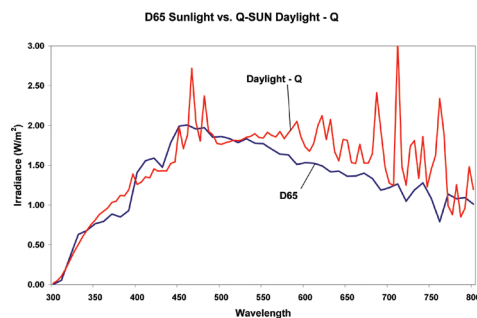

This chart shows the Q-SUN chamber's Daylight-Q filter compared to D65 outdoor sunlight.

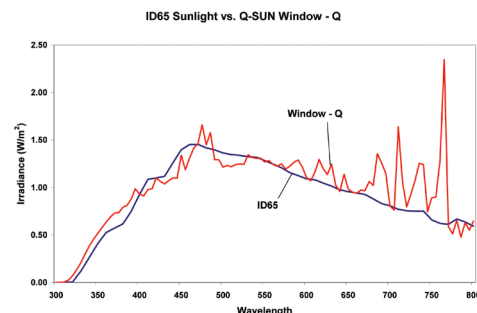

This chart shows the Q-SUN chamber's Window-Q filter compared to ID65 indoor daylight.

Two Sizes - All the Critical Features

The Q-SUN tester comes in two sizes with all the critical features for pharmaceutical and cosmetics applications. The Xe-1 is available in a convenient tabletop size and offers an optional chiller. For labs that need high throughput, the Xe-3 is a free-standing unit with triple the capacity. The Xe-3 also provides precise control of relative humidity.

ICH Guidelines

Photostability testing of drug substances and products can be done several ways according to the ICH Guidelines. Regardless of whether Option 1 or Option 2 is chosen or which light source within each option is used, the exposure requirements are defined as **minimum** values within the visible and ultraviolet ranges:

- Visible Light: 1.2 million lux-hours.
- Ultraviolet Light: 200 w-hr/m² of integrated UV energy (spectral range of 300-400 nm)

Under the ICH Guidelines, it is perfectly acceptable to expose specimens to higher dosages in either range. In fact, light sources that meet the D65 or ID65 definition will achieve the visible light energy dosage well before the ultraviolet energy dosage. Therefore, a D65 or ID65 light source will create an “overexposure” of uv energy upon completion of the test.

Using the Q-SUN Test Chamber to Meet ICH Guidelines. The Q-SUN tester can be an ideal solution for ICH Guidelines Option 1 because of its excellent simulation of sunlight. Because the light intensity is not specified in ICH, a variety of test settings can be run in the Q-SUN unit to meet the minimum exposure requirements.

The ICH Guidelines do not specify whether the spectrum should match D65 or ID65. Commonly, the ID65 spectrum is simulated using a Window-Q filter in the Q-SUN because it represents sunlight entering a building interior. However, a daylight filter is appropriate for outdoor applications such as hair dyes and sunscreens.

The table below gives some examples of Q-SUN chamber set points using a Window-Q filter and a 420 nm control point. Other control points and filters may be used - contact Q-Lab for details.

Set Point @420 nm	UV Irradiance (300-400nm)	Illuminance	Test Duration	UV Dosage	Visible Light Dosage
0.62 W/m ²	28.3 W/m ²	50.4 klx	24 hours	673 W-hr/m ²	1.21 million lux-hr
1.10 W/m ² *	49.7 W/m ²	89.5 klx	13.5 hours	671 W-hr/m ²	1.21 million lux-hr
1.47 W/m ²	66.5 W/m ²	119.6 klx	10 hours	665 W-hr/m ²	1.20 million lux-hr

* 1.10 W/m² is a common set point in photostability tests.

Temperature and Relative Humidity Control

ICH Guidelines do not include temperature or relative humidity requirements. Tests are commonly performed in a variety of temperature environments, including normal room temperatures. Adding the optional chiller to the Q-SUN unit is recommended to achieve common room temperature conditions in the test chamber.

Although chamber relative humidity may not be critical for testing drug substances enclosed in packaging or test bottles, the Q-SUN Xe-3 provides precise control of RH for tests that may benefit.



Bottles and other oddly-shaped samples can be easily tested in a Q-SUN chamber due to its flat specimen tray and variety of mounting options.

Q-Lab Corporation _____ www.q-lab.com



Q-Lab Headquarters
Westlake, OH USA
Tel: +1-440-835-8700
info@q-lab.com

Q-Lab Florida
Homestead, FL USA
Tel: +1-305-245-5600
q-lab@q-lab.com

Q-Lab Europe, Ltd.
Bolton, England
Tel: +44-1204-861616
info.eu@q-lab.com

Q-Lab Arizona
Buckeye, AZ USA
Tel: +1-623-386-5140
q-lab@q-lab.com

Q-Lab Deutschland GmbH
Saarbrücken, Germany
Tel: +49-681-857470
vertrieb@q-lab.com

Q-Lab China 中国代表处
Shanghai, China 中国上海
电话: +86-21-5879-7970
info.cn@q-lab.com